

Duplex Vertical Tank Mount

1.5 HP, 2HP & 3 HP



Description

NFPA 99 compliant duplex lubricated rotary vane medical vacuum package. One piece vertical tank mount system is designed for single point connections. All systems are thoroughly tested prior to shipping. This system ships with:

- Inlet and discharge threaded flex connectors
- Hard copy and CD of O&M manuals

Package Details

Duplex medical vacuum system, consisting of two continuous duty, lubricated rotary vane, air cooled vacuum pumps. The vacuum pumps, vacuum receiver and control panel shall be factory assembled in a one piece package, pre-plumbed and pre-wired for single point connections. The package includes: 5.0 micron inlet filters, vibration isolation feet under each pump, exhaust and intake manifold plumbing with flexible metal vibration isolators, pump isolation valves, 3 valve by-pass and main system shutoff valve with flexible metal isolator for connecting to piping system. The package meets or exceeds NFPA 99 current edition.

Lubricated Rotary Vane Vacuum Pump

Each medical vacuum pump shall be single stage, rotary vane, continuous duty type with a recirculating oil system filtered by an easily serviced spin on oil filter. Each pump shall be capable of operating continuously through its full range up to 29.3" HG. Each pump shall contain an internal mist eliminator providing 99.9% oil free exhaust and a drip leg with drain valve on pump exhaust discharge. Each vacuum pump is equipped with an inlet check valve and anti-suck back valve, vacuum inlet filter with a 5.0 micron replaceable element and a temperature safety switch shall be installed at the exhaust port of each vacuum pump. The vacuum pumps shall be filled with a high quality extended life, full synthetic lubricant. Each vacuum pump is air cooled and no water shall be required for normal operations.

Electric Motor/Drive

Each vacuum pump shall be direct driven by a 1,750 RPM C-face, totally enclosed, fan cooled, TEFC electric motor. Motor shall be 208-230/460 volt, three phase with a service factor of 1.15.

Vibration Isolation

Each pump and motor assembly shall be mounted to the base frame/tank with seismically restrained vibration isolators.


Standard Control Panel:

The system includes a duplex control panel in a NEMA 12 enclosure. The control panel is UL 508 listed and is designed to operate an NFPA 99 compliant medical vacuum system. The panel includes: individual through-the-door disconnects for each motor, motor circuit breakers, magnetic starters, adjustable motor overload protection, dual control voltage transformers, alternating circuit, hand-off-auto switches, auto/lead pump selector switch, power on light, high temperature lights, transformer failure light, pump run lights, vacuum switches and alarm contact connection points. The control panel will start up immediately after power failure and no boot-up time is required.

Optional

- Galvanized tank.
- Water cooled oil coolers for high ambient conditions.
- Touch-Screen Interface Control System (HMI/PLC) with built-in ethernet, system monitoring, maintenance alerts and multiple protocols for facility automation systems.

Lubricated Rotary Vane Vacuum Package Duplex Tank Mount 1.5, 2 HP & 3 HP

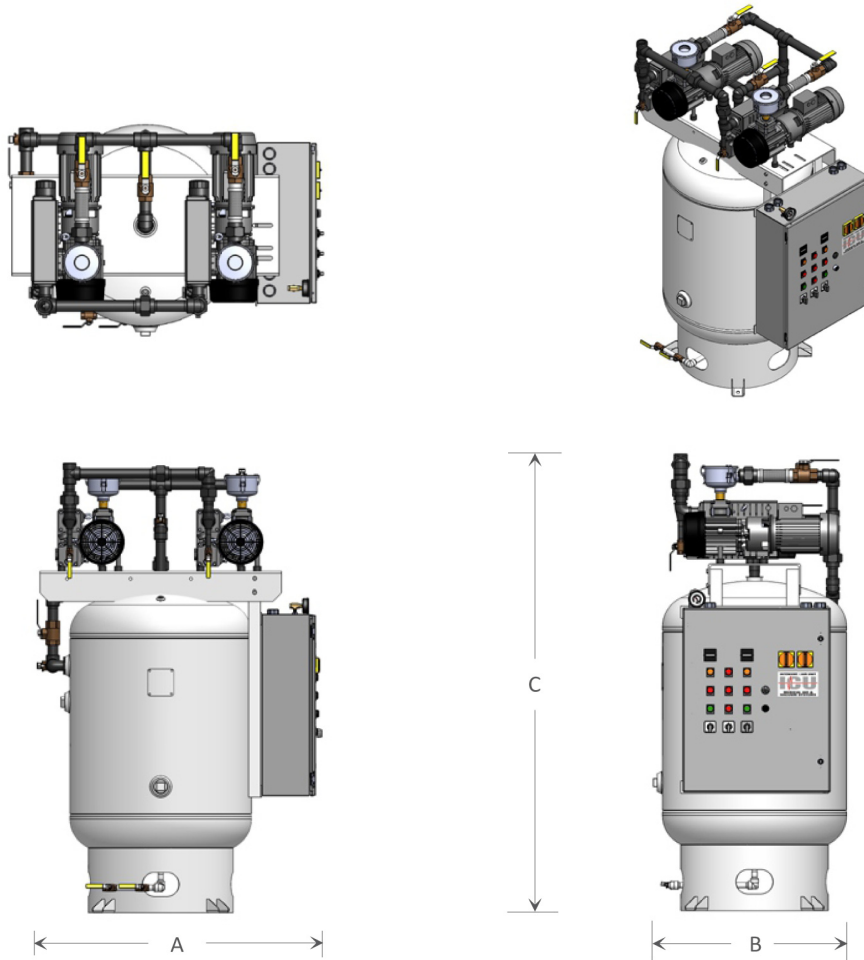

There are recommended minimum clearances for equipment maintenance and ventilation.
Please contact us toll free at 800.333.0811

Shown: Duplex 2 HP, 120 gallon, vertical tank mount, standard lubricated rotary vane configuration.

Model Number	Tank Gals	Horsepower		SCFM @ 19" HG ^(1,2)		dB(A) ⁽³⁾ Level	System FLA			BTU HR	Dimensions			NPT In	NPT Out	Weight (LBS)
		each	pkg	each	pkg		208V	230V	460V		A	B	C			
V73M-DTB-3VS	80	1.5	3.0	7.0	7.0	70	9.0	8.4	4.2	3,817	42.0	28.0	67.0	1.25	1.25	706
V73M-DTB-4VS	120	1.5	3.0	7.0	7.0	70	9.0	8.4	4.2	3,817	47.0	35.0	74.0	1.25	1.25	922
V73M-DTC-3VS	80	2.0	4.0	7.0	11.0	70	11.2	11.2	5.4	5,089	42.0	28.0	67.0	1.25	1.25	734
V73M-DTC-4VS	120	2.0	4.0	11.0	11.0	70	11.2	11.2	5.4	5,089	47.0	35.0	74.0	1.25	1.25	950
V73M-DTD-3VS	80	3.0	6.0	17.0	17.0	70	16.8	15.6	7.8	7,633	42.0	28.0	67.0	1.25	1.25	838
V73M-DTD-4VS	120	3.0	6.0	17.0	17.0	70	16.8	15.6	7.8	7,633	47.0	35.0	74.0	1.25	1.25	1,034

1. Vacuum pump data is shown with one pump in reserve. 2. Pump capacities are shown in SCFM (Standard Cubic Feet per Minute).
3. dB(A) levels are shown as only one pump running. 4. Data is shown with equipment running at normal operating conditions with a maximum ambient temperature of 105°.

